
ELECTRIC HYDRAULIC HIGH SHEAR HOMOGENIZER

For Processing Samples as Low as 3 mL

New and Improved HMI!

- Increased flow rate
- Superior cooling through specially designed product heat exchanger



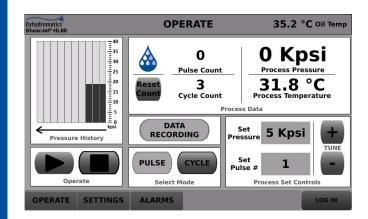

The ShearJet HL60 electric hydraulic processor is a unique R&D model designed to process samples as small as 3 mLs in distinct passes or continuous runs. Multiple samples may be processed in a short time span. Small dead volume (due to fewer high pressure connections) and self-priming leads to reduced waste. The HMI touchscreen display puts simplicity at your fingertips, displaying relevant information in real time.

Performance:


- Multiplier pump creates unique pressure profile
- Specially designed product heat exchanger ensures clean, superior cooling
- Proprietary fixed-geometry diamond Reaction Chamber[®] module with specially designed micro channels achieves superior, consistent, uniform, repeatable mean particle size reduction
- Durable, robust design allows for continuous operation
- Data acquisition via USB or Ethernet

Easy-to-clean:

Minimum connections allow for easy cleaning

Economy:

- Save on expensive APIs
- Process samples in minutes
- Low maintenance cost
- Minimal seals and metallic connections
- Diamond chamber for longer life

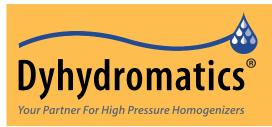

Q904-0021-R012 5/10/24 All descriptions and technical data subject to change.

www.dyhydromatics.com

Specifications

Minimum sample size:	3 mL (minimum with cooling jacket 5 mL)	
Working sample size:	3-500 mL and up	
Pressure:	5-30 Kpsi (345-2068 bar)	
Stroke:	1 mL	
Output*:	64 mL/min @ 30 Kpsi (2068 bar)	
HMI:	7″Touch screen	
Power:	Electric hydraulic 12A @ 220V	
	Ask about 110V model	
Dimensions:	28L" x 16W" x 24H" (71,1 x 40,6 x 61,0 cm). Small foot print saves valuable laboratory space!	
Noise level:	68 decibels; more quiet than competitive equipment	
Operating time:	System can run anywhere from seconds to hours	
	*Dependent on Reaction Chamber module size.	

*Dependent on Reaction Chamber module size. Buna-N is a registered trademark of E. I. DuPont. Luer-Lok" and "Luer-slip" are registered trademarks of Becton Dickinson.

Features


Features for HL60	Standard Model	Premium Model
Glass reservoir 100 mL with sanitary connection	۵	
Stainless steel reservoir 100 mL		۵
Inlet thermocouple		۵
Two 10 mL syringes	۵	۵
Luer lock connections	۵	۵
Ceramic plunger	۵	۵
Inline pressure transducer		۵
Diamond Reaction Chamber® module for longer life	۵	۵
and wear resistant to abrasive applications		
Jacketed helical heat exchanger	۵	۵
Touchscreen HMI with advanced features	۵	۵
Tools & spares	۵	۵
Chiller for cooling product heat exchanger	option	option
Gaskets and O-rings in TFE, Kalrez [®] , Buna-N [®] , EPDM, or EPR	option	option
Pharma-grade product contact parts	option	option
FAT/SAT	option	۵
IQ/OQ	option	option
Onsite installation and operator training	۵	۵

HMI

- · Easier to use touchscreen
- Improved Data logging
- Export data recording via USB or FTP Client
- Access levels
- Historic data graph
- Calibration input
- Allows operator to set process duration in pulse mode. System stops automatically when completed.





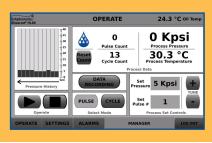

Dyhydromatics Corporate Headquarters

63 Great Road, Suite 103 Maynard, MA 01754 Phone +1 978 461 0200 • Fax: +1 978 461 0201 Customer Service: (Option 4) info@dyhydromatics.com Sales: sales@dyhydromatics.com

Dyhydromatics European Head Office

Kelterstrasse 9, 75236 Kaempfelbach, Germany Customer Service Europe: Phone +49 7231 167 4048 Sales: info-europe@dyhydromatics.com

Dyhydromatics, India

A-7, Mangal Society, Shreenagar Wagle Estate, Thane (W) - 400604, India Phone +91 77381 42553 Email: info-india@dyhydromatics.com FORCE BLUE is the only nonprofit organization in the world that provides "mission therapy" for former combat divers by retraining, retooling and deploying them



ooling and deploying them on missions of environmental conservation, preservation and restoration. We are pleased to support FORCE BLUE by donating a portion of proceeds from every new equipment sale.